



Caution: Federal law restricts this device to sale by or on the order of the physician.

The Personal KinetiGraph (PKG) is intended to quantify kinematics of movement disorder symptoms in conditions such as Parkinson's disease, including tremor, bradykinesia and dyskinesia. It includes a medication reminder, an event marker and is intended to monitor activity associated with movement. The device is indicated for use in individuals 46 to 83 years of age.

PKG® is Personal KinetiGraph® in the USA.
PKG® and KinetiGraph® are registered trademarks of Global Kinetics Pty. Ltd.

Manufactured by:
GKC Manufacturing Pty Ltd
Level 9, 31 Queen Street, Melbourne, VIC, 3000, Australia


Distributed by:
Global Kinetics USA Corporation
10913 Bren Road East, Hopkins, MN 55343

info@globalkineticscorp.com | www.globalkineticscorp.com

MM0154-US Rev 1.0 March 2020



Wearing your
PKG® Watch

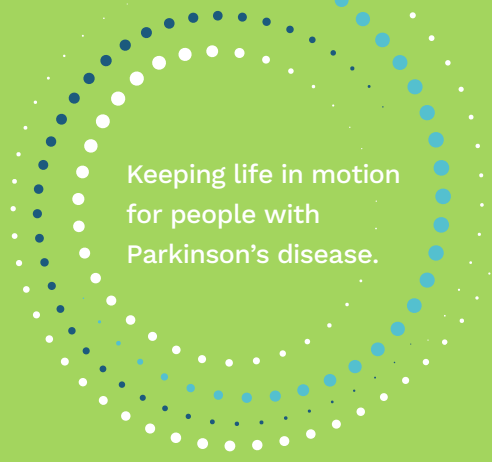

Keeping life in motion
for people with
Parkinson's disease.

Your doctor has referred you for a Personal KinetiGraph® (PKG)

The PKG Watch is a wrist worn device that records your Parkinson's symptoms and reminds you when to take your Parkinson's disease medication as prescribed by your doctor.

This recording is done at home over 6–10 days while you go about your day-to-day activities.

The information collected by the PKG Watch can tell your doctor more about your movement during the day, including if you are moving a bit slowly or having difficulty in performing movements. The PKG can help your doctor assess the current state of your Parkinson's disease and assist them in making a decision about the right type and amount of medication for you at the right time.



Watch

a video on using the PKG Watch at bit.ly/thepkgus



Educate

yourself with the instructions they will send you.



Act

now and wear the PKG Watch as soon as you receive it.



Return

the PKG Watch promptly so we can send your results to your doctor.



Your PKG Road Map

Your PKG journey is simple and easy to navigate. Follow your doctors instructions and make sure to wear the PKG Watch as soon as it arrives.

Order Placed

Your clinician has ordered you a PKG Watch.

____/____/____
Date Ordered



Prepare

Watch the video on "How to Wear the PKG Watch," go to, bit.ly/thepkgus

Customer Service Call

A representative will call to confirm your medications. Make sure to answer the phone!



Your Watch Arrives

Wear your PKG Watch immediately for the next 7 days. Your provider is waiting on your data



Return Your Watch

The process is almost complete! After the 7th day, use the pre-paid shipping label and drop in any mailbox.



Your Data

Global Kinetics specialists are analyzing your data!



Congratulations! You have reached your destination!

Your clinic receives your PKG report.



Get more information about your symptoms to optimize your Parkinson's treatment.



FAQ's

What is Objective Measurement?

Your clinician uses many tools to help you manage your Parkinson's symptoms. Objective measurement, using the PKG, quantifies your movement symptoms over the course of a week while you live your life. This information is especially useful as it relates to your current medications which is why we confirm with you and provide medication reminders on your PKG Watch.

How will my doctor use my PKG data?

Your PKG shows your doctor how your movement changes though the course of the week and your movement symptoms response to your medications. Along with a detailed clinical assessment and conversation with you, your doctor may chose to make adjustments to your medications, maintain your current treatment, or suggest alternative therapies. Your PKG is not a substitute for expert clinical care.

How common is PKG for Parkinson's?

While PKG may be new to you, over 50,000 PKGs have been performed to date, in the United States, Europe and Australia.

Can my PKG Watch get wet?

Your PKG watch is designed to be water resistant. You can wear it in the shower, but remove it before activities like swimming or bathing.

What do I do if I have questions?

Your clinic is the best place to call if you have any questions related to your medical care. If you have questions related to your PKG Watch, please visit our website, email or call. We're here for you!

Contact PKG Patient Support with questions.

Email: pkgHelpUSA@GlobalKineticsCorp.com

Phone: 952-873-1737, 8:00 a.m. to 6:00 p.m. CST